

■ ■ wissner-
■ ■ bosserhoff

Functional elegance in nursing care



sentida

Universal low nursing beds





Safety and comfort in perfection

sentida – universal low nursing beds

Falling is one of the greatest risks for the physical health and psychological well-being of elderly people. For this reason, systematic fall prevention is of crucial importance in care for the elderly, making it an important quality feature of a nursing home. The consequences of falls are often very serious for elderly people and, in the worst case, even irreversible.

The fear of further falls is one of the most common psychological consequences, leading to limited physical activity, a loss of self-confidence and social isolation. With the sentida range of low nursing beds, wissner-bosserhoff has achieved a consequent and above all uncompromising implementation of the requirements of modern nursing care: Maximum safety hand in hand with a cozy design.

Innovative solutions such as the 3-stop strategy and the patented SafeFree® side guard concept reduce falls as well as the consequences associated with them and set standards in terms of safety, mobility, height adjustability and coziness.

5

million
falls of elderly
people per year:

Almost one third of people aged 65 and above, as well as 50% of those aged 80 and above fall at least once a year.

In Germany alone, there are about 5 million cases of elderly people falling each year.

(Source: "Gesundheit und Krankheit im Alter" (= health and illness in old age).
Publisher: Karin Böhm, Federal Statistical Office of Germany / Clemens Tesch-Römer, German Centre of Gerontology / Thomas Ziese, Robert Koch Institute, Berlin 2009)



Prevention
with sentida
low nursing beds

for use
in institutions with
high standards of universal
care close to the ground:
Nursing homes and
homes for the elderly,
care resorts, palliative wards,
dementia and geriatric wards
and rehabilitation
institutions.

Living concept

sentida | Design G10



The selection of materials and colors follows current, natural living trends. This way, a cozy, warm atmosphere, in which the resident feels safe and protected, is created.

sentida Portfolio

The entry model



The mobile one



The flexible one



The cozy one



The extra mobile one



Model	sentida 1	sentida 3	sentida 4	sentida 5	sentida 6
Movability	In low position	At every bed height	At every bed height	At every bed height	At every bed height
Castors	Two fixed castors at the foot end Two 100 mm castors at the head end	Four 50 or 75 mm castors at the foot end double castors Two 100 or 160 mm castors at the head end fixed castors	50 or 75 mm double castors	50 mm double castors, completely covered by cozy undercarriage cladding	100, 125 or 150 mm easy-roll castors or 125 mm double or design castors
Braking capacity	Electric brake with handset	Electric braking with handset foot-end brake bar	Head-end and foot-end brake bar	One pedal brake for each axle with 2-castor central braking system	Foot-end brake bar with 4-castor central braking system
Safe working load	225 kg	225 kg	225 kg		
Height adjustment	25-80 cm	25-80 cm	25-80 cm	27-80 cm	27-80 cm with 100 mm easy-roll castors 29,5-82,5 cm with 125 mm easy-roll castors 28-81 cm with 125 mm double or design castors 32-85 cm with 150 mm easy-roll castors



Benefits at a glance



Protection and mobility without FDM*

SafeFree® – the patented side guard concept

| 12 - 13

Safety by day and night

SafeSense® – the intelligent bed-exit sensor system

| 14 - 15

Safe nursing care, mobilization and sleep

3-stop strategy

| 16

Flexible and time-saving

3 bed lengths that can be adjusted without tools

| 17

Effective bedsore prevention

Double retraction | 3 patient surfaces to choose from

| 18 - 19

Intuitive operation

Clearly structured and ergonomic

| 20 - 21

Energy efficiency and stability

Thanks to innovative drive concept

| 22 - 23

We think ahead

memoriana | Room concept for dementia patients

| 24 - 25

Additional options

Design | Colors and decors | Features and accessories

| 26

Living concept

sentida | Design Ts10



The selection of materials and colors follows current, natural living trends. This way, a cozy, warm atmosphere, in which the resident feels safe and protected, is created.

Protection and mobility without FDM*

SafeFree® the patented side guard concept

The SafeFree® divided side guard system is patented and voluntarily certified. The technology has been tested in practice more than 500,000 times. You can lower the space-saving side guards quickly in only 2 seconds and individually adjust them to four heights, depending on the care situation and size of the resident.

1
No protection needed



Lowered, the divided side guards look homely and provide unhindered access to the resident for daily routine nursing tasks. No need to take off and store accessories that are in the way, such as a middle post or fixed side guards.

2
Soft protection (no FDM*)



For safety at night, it is often sufficient to combine the low position of the bed with the first height level of the head-end side guard. This avoids unnecessary barriers and effectively reduces the risk of falls on the long term. The resident can orient himself without being restricted.

3
¾ protection (no FDM*)



Creates a deliberate bed egress gap at the foot end to reduce freedom-depriving measures (FDM) while at the same time maintaining a high safety standard.

4
Full protection (FDM*) for mattresses up to 18 cm




The highest position provides comprehensive protection of the resident. The distance between the two parts of the side guards is minimal, making them just as safe as a continuous side guard and exceeding the requirements of the standard IEC 60601-2-52. It is not necessary to stock up on accessories to close a central gap when using an abdominal belt for restraint.



If you use a thicker mattress, you can raise the side guards by an additional 6 cm.


Safe mobilization

Divided side guards are very useful for effective mobilization, as they are easy to grab and clasp. They can be used as a safe standing-up aid and verifiably support the resident's mobility.

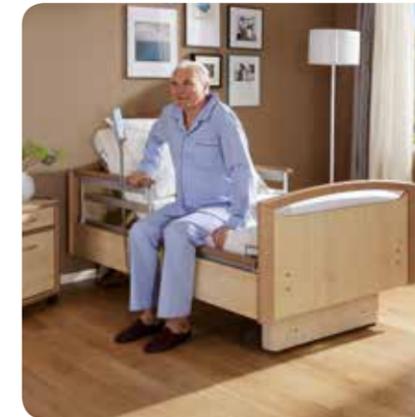


All-round access for nursing tasks = with 4 divided side guards

- Flexible use of the bed for standard care, even in limited spaces
- Flexible access from both sides of the bed in case of increased need for care
- Flexible use of the bed for bed egress on the left or right side to mobilize the resident
- Avoiding follow-up costs for upgrading divided side guards

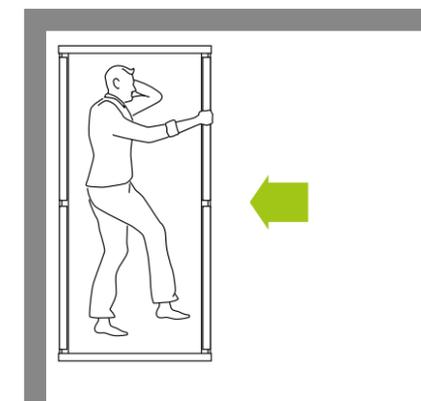

Standard nursing - access from one side is sufficient

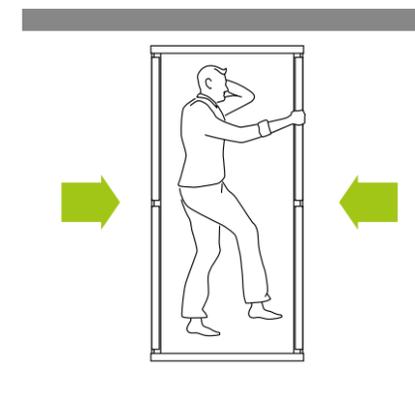

Increased need for care - access from both sides needed

More safety due to smaller central gap* < 60 mm

- Targeted compliance with the standard IEC 60601-2-52
- Safety in working position, i.e. also with raised backrest or leg rest
- Safety without using an additional central gap protector as accessory for restraint with an abdominal belt in accordance with the BfArM** (German Federal Institute for Drugs and Medical Devices)

tested in practice

500.000 +

times

Videos on the side guard concept

Safety by day and night

SafeSense® – the intelligent bed exit sensor system

wissner-bosserhoff, your reliable partner for the reduction of freedom-depriving measures (FDM):
With SafeSense®, wissner-bosserhoff presents an intelligent assistance sensor system, which reliably helps reducing freedom-depriving measures day and night. It is equipped with an automatic night light, which is activated when the resident leaves the bed and is deactivated when the resident is back in bed again. SafeSense® is compatible with standard nurse call systems and can be retrofitted to wissner-bosserhoff nursing home beds manufactured in 2002 or later



↓

More safety –
less freedom-
depriving measures
(FDM)

- Faster response time in the event of falls
- Greater safety during night and weekend shifts
- Ideal in combination with sentida low nursing beds and SafeFree® side guards

Automatic night light

- Night light is automatically switched on when the resident leaves the bed
- Night light is switched off again as soon as the resident is back in bed

Individual alarm system


- Freely adjustable immediate alarm from 0 seconds to 30 minutes
- Adaptation to the resident's mobility



Video SafeSense®

Integrated or available as retrofit kit

- Accessory option for current WIBO nursing home beds
- Retrofit kit for electrical WIBO nursing home beds manufactured in 2002 or later
- Available with cable or wireless

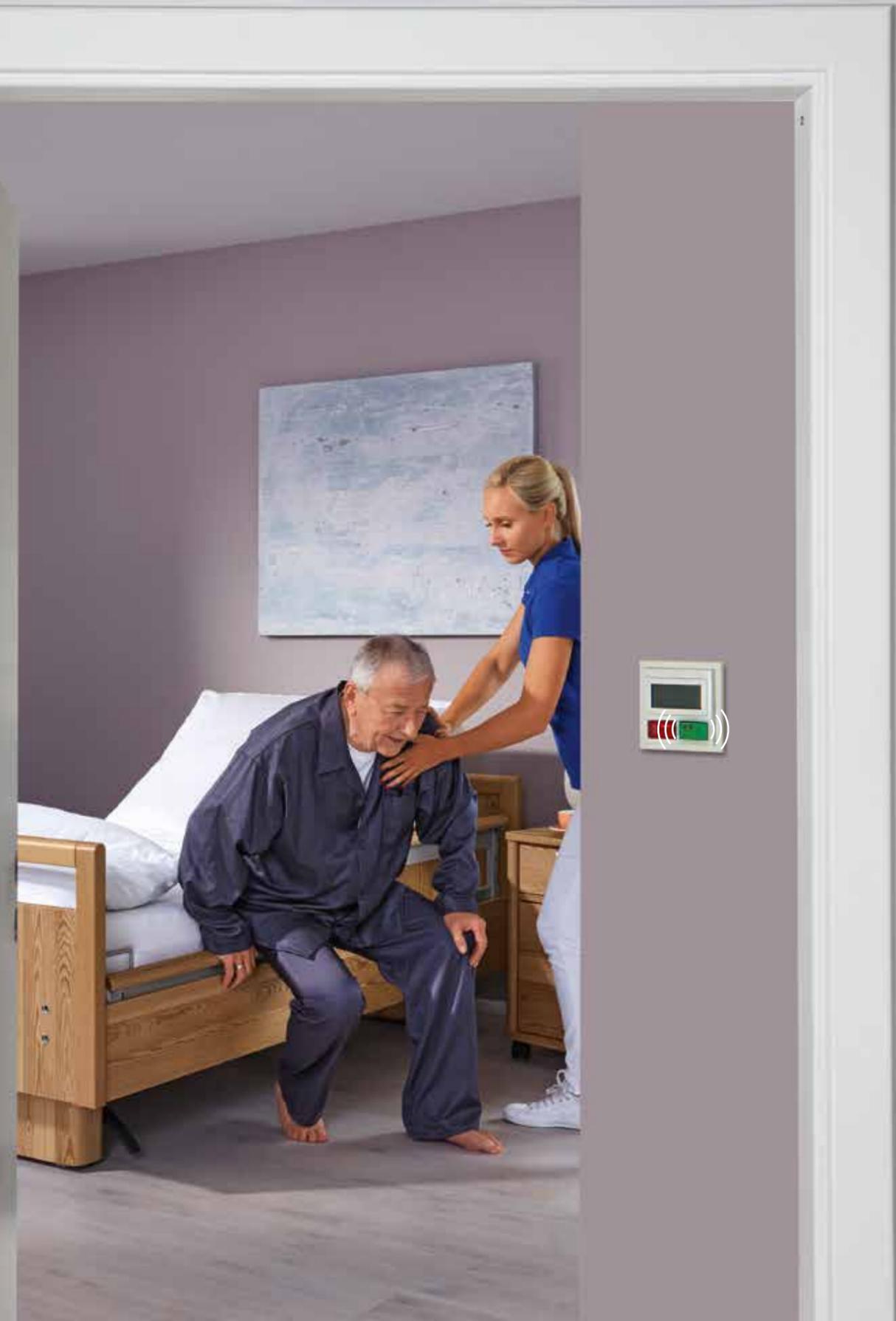

Safe nursing care, mobilization and sleep

3-stop strategy

The **IVS intelligent verticalization system** not only features a great **height adjustment range of 27 to 80 cm**, but due to its integrated control sensors also prevents collisions when setting the bed to the comfort seating position. This becomes an indispensable advantage, especially when combined with the bed extension to 220 cm.

80 cm

Optimal nursing position




41 cm*

Safe mobilization at chair height



27 cm*

Fall prevention during sleeping phases

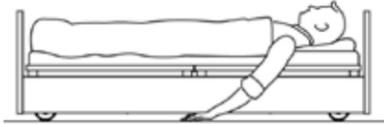



IVS – intelligent verticalization system

- Comfortable adjustment range from 27 to 80 cm
- First automatic stop at 41 cm at chair height for bed egress
- Collision protection when adjusting the bed from the low position and when using a bed extension


Flexible and time-saving

3 bed lengths that can be adjusted without tools

All standard wissner-bosserhoff nursing home beds are equipped with an integrated bed extension. With the **express bed extension** (see image below/option), sentida - **with divided side guards** - can be adjusted flexibly and without tools to different bed lengths. The requirements of the standard IEC 60601-2-52 are complied with in each configuration.


200 cm

Bed length of 200 cm:
The divided side guards offer protection without a central gap > 318 mm* and without a protector

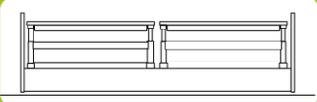


210 cm


Alternative bed length of 210 cm:
Protection with longer side guards at the foot end without accessories and without a gap > 318 mm**

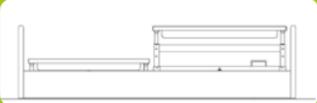

220 cm

Bed extension to 220 cm:
is temporarily adjustable with divided side guards for sentida - complete protection through foot-end protector (accessory) is guaranteed

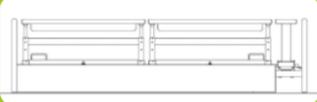

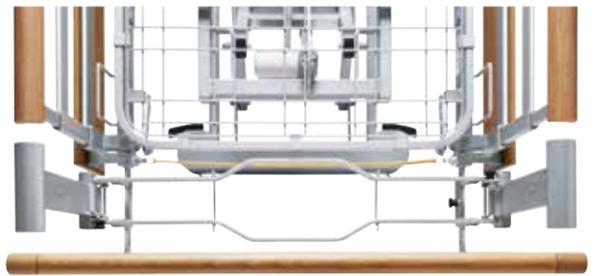

Bed extension with express release (option)

Effective bedsore prevention


Double retraction | 3 patient surfaces to choose from

Patient surface ergonomics


The ergonomic division of the patient surface helps prevent pressure and therefore contributes to pressure ulcer reduction. The dimensions of the 4-part patient surface of wissner-bosserhoff nursing home beds are based on findings from anthropometrics, i.e. on the average size and height of a person.

On the basis of professor Eigler's findings, it can be seen that the backrest area should be at least 80 cm long in order to effectively support the cervical spine. An additional reduction of pressure can be reached with the BiNetic patient surface for patients with above-average height due to its

integrated bed extension. The patient surface can be extended to 220 cm either temporarily and without tools using the express bed extension, or permanently with the bed extension set, preventing the legs from bumping on the foot end.

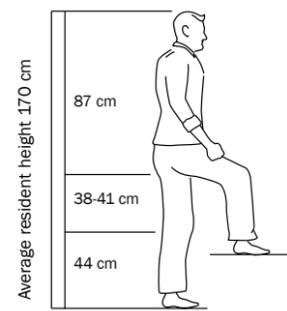

Anthropometry according to Prof. Eigler *

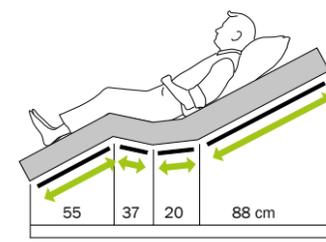


BiNetic patient surface 200 cm (long backrest)

Patient surface modularity


You can choose between 3 different patient surface options:

- 01 The wire mesh patient surface Aero ensures good aeration.
- 02 The synthetic mats of the Easyclean patient surface are especially easy to clean.
- 03 The comfort patient surface with its 55 spring elements not only promotes bedding comfort but also makes a valuable contribution to pressure ulcer prevention.

In general, the patient surface partitions can be easily removed and cleaned.

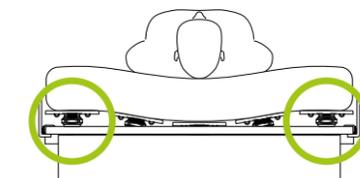

Comfort patient surface with 55 spring elements

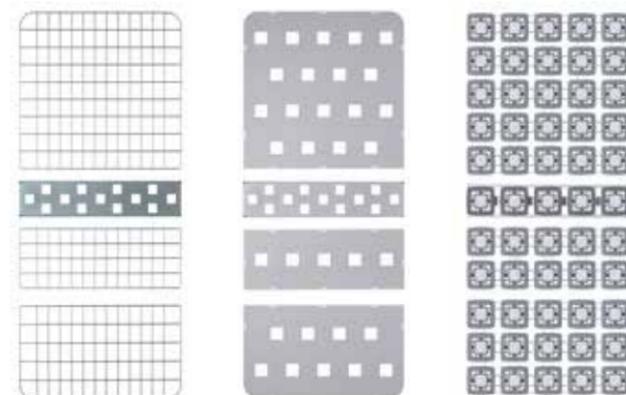

01

02

03



Detail comfort patient surface


Effectiveness of double retraction confirmed*

For years now, wissner-bosserhoff has been pursuing a practicable concept for reducing the risk of pressure ulcers: The double retraction as a component of the BiNetic patient surface. With the double retraction, not only the backrest but also the leg rest is located further apart from the seat section, so that

the surface is enlarged. Today, results from pressure measurements (thermal imaging) clearly show that there is less pressure with double retraction than with single retraction and significantly less pressure than with the standard patient surface.




WIBO patient surface with notable relief

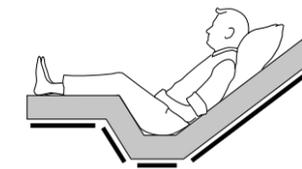

6 cm ← → 11 cm

Double retraction = pressure ulcer prevention

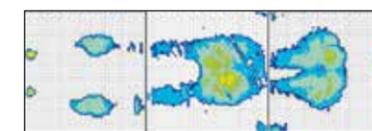


Thermal image of a BiNetic patient surface



Standard patient surface without relief



Single retraction = risk of pressure ulcer

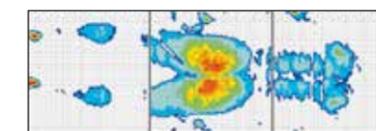

Thermal image of a standard patient surface

Up to **2.7** million people in Germany have chronic wounds such as bedsores. Not least due to the demographic change, elderly people suffer from chronic wounds particularly often. The risk for people over 90 is at more than 11%. Pressure ulcers are very unpleasant and straining for residents. Therapies usually take several weeks or months and may often cause further discomfort.

Source: BVMed 2012 study


Double retraction = **17cm** more space in bed

Intuitive operation


Even more straightforward and ergonomic

Already at first glance, the handset of the sentida nursing home beds impresses with its clarity and ergonomic form: fewer buttons and large symbols ensure that operation can be learned intuitively.

The pre-selection function makes the handset single fault safe – triggering a function by accidentally pressing a key is not possible.

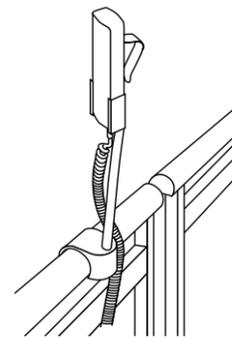

The handset is easy to reach for both the resident and staff at all times (accessory handset holder, model No. 10-0419 and 10-0421)

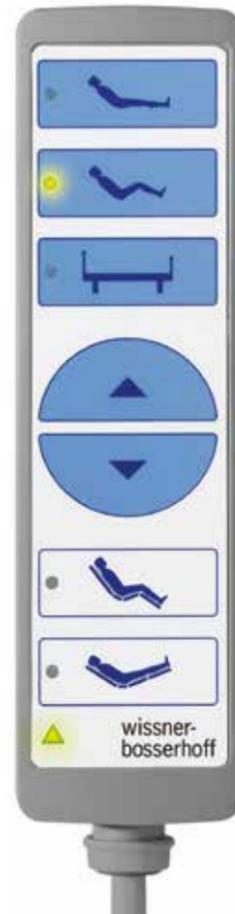

Backrest adjustment

Auto contour


Height adjustment

Up

Down

Automatic comfort seating

Automatic positioning

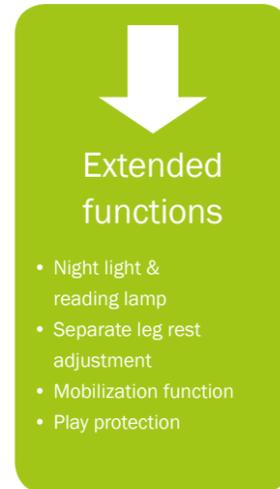

- Night light & reading lamp
- Separate leg rest adjustment
- Mobilization function
- Play protection



Video tutorial, also at www.wi-bo.com/de/Ueberuns/Anleitungsvideos

The automatic positions of sentida, adjustable at the push of a button, increase comfort when lying in bed and support pressure ulcer prevention for the benefit of the resident while reducing the strain on nursing staff.

Watching television, reading or eating breakfast in bed with automatic comfort seating

The automatic comfort seating brings the bed into a pre-selected seating position for watching TV and reading comfortably, for eating or communicating at eye level. Mobilization using the lifter is effectively supported in this position. In addition, this position facilitates breathing and digestion. It trains spatial orientation and stimulates circulation for residents confined to bed. The double retraction provides additional relief of pressure in the sacral area.



3 safety levels

The locking function is integrated into the back of the handset, meaning that you do not have to look for a control box at the patient surface panel or bend down.

The three handset modes can be selected with a non-magnetic key:

- Care mode:** All selection buttons are unlocked.
- Locking mode:** All functions of the handset are locked.
- Resident mode:** The resident can use the upper three buttons.



Key for safety levels (magnet-free)

Bedsore prevention and keeping heels free with automatic positioning

The preset legs up position promotes decongestion with venous diseases and at the same time supports pressure ulcer prevention. The Semi-Fowler 30° position with knee lifting guarantees a reduced decubitus risk* thanks to a favorable weight distribution in the sacral region, which is especially at risk of developing pressure ulcers.



Energy efficiency and stability

Thanks to innovative drive concept

With the new motor, the sentida 5 and 6 nursing home beds are even more user-friendly, stable and energy-efficient. The integrated LED technology, energy-saving features and low-maintenance design are vital requirements for the cost-effective operation of a nursing home.



Even more user-friendly

The new sentida 5|6 is equipped with an LED bedside light connection integrated in the motor. You only need one single cable from the bed to the power supply. The LED light can also be switched on and off using the handset. You can now also watch video tutorials on how to operate the bed by scanning the QR code on the bed.



Video tutorial, also at: www.wi-bo.com/de/Ueberuns/Anleitungsvideos

Stable and safe

The new motor is especially powerful. Its performance is especially impressive under high load. It has a safe working load of up to 270 kg (sentida 5|6) and is suitable for a maximum resident weight of up to 235/255 kg (depending on the accessories used).



Economical and energy-efficient

With continuously increasing electricity prices, energy efficiency is an increasingly important factor. With an extremely low stand-by consumption of only 0.5 W, the new motor is especially economical. A clever cable management reduces wear and tear. The power cable can be exchanged separately, causing less storage costs.

The following calculation example demonstrates how the stand-by consumption of nursing home beds can influence the costs:



Product	Stand-by operation	Electricity costs over the course of 15 years*
The NEW sentida	0.5 W	€1,950
Bed A	3.76 W	€15,345
Bed B	2.5 W	€10,200
Bed C	1.7 W	€6,900

*Example: Institution with 100 beds

soluna

Designer reading light

The elegant shape of the partially translucent lampshade stands out as a special design element of soluna lights. Equipped with a swivel arm, you can precisely direct the light towards the object to be illuminated and it blends in harmoniously with the nursing home beds' ergonomics. Benefits of using LEDs: a long life, low energy consumption, no heat and absolutely safe operation.



With its variable light head and the individually adjustable light focus, soluna is ideal both as a reading light and bedside light.



soluna with undulating lampshade



soluna with smooth lampshade



We think ahead

memoriana | Room concept for dementia patients

Adequately taking care of people with dementia is a particular challenge for nursing staff and relatives. In cooperation with dementia experts from Great Britain and Germany, wissner-bosserhoff has developed a room concept especially for dementia patients. Specifically adapted to the needs of the patient, the new concept offers a considerable reduction in the strain on those providing care. The focus of the development includes three targets for the care of residents with dementia:



Improving orientation

The resident is offered reference points for orientation in time and space, as well as colored orientation aids. The targeted integration of colors, forms and objects helps in terms of spatial orientation – worries and frustrations are reduced.



Increasing independence

The more independently the resident can move in his surroundings, the less strain on the nursing staff. Functional furnishing elements such as handrails and open storage areas promote the dementia patients' mobility and independence. Other areas of the room are, however, only accessible to those providing care.



Maximizing safety


Resident safety is of utmost importance. Low nursing beds with side guards, handrails on furniture and walls as well as visual orientation aids provide the person diagnosed with dementia a maximum level of safety without feeling restrained.

100

different
causes

Dementia can have up to 100 different causes. As complex as the causes and triggers of dementia is the course of the illness, which may differ from resident to resident.


Source: www.bmg.bund.de
German Federal Ministry for Health

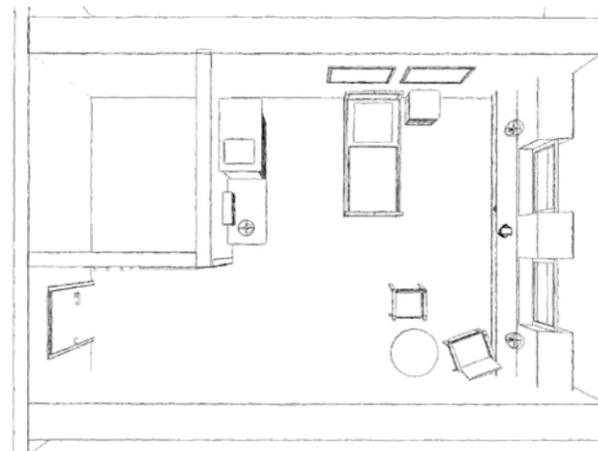

For more detailed information on the dementia room concept, have a look at the separate memoriana brochure or contact us by e-mail: info@wi-bo.de




Design | head / foot sections

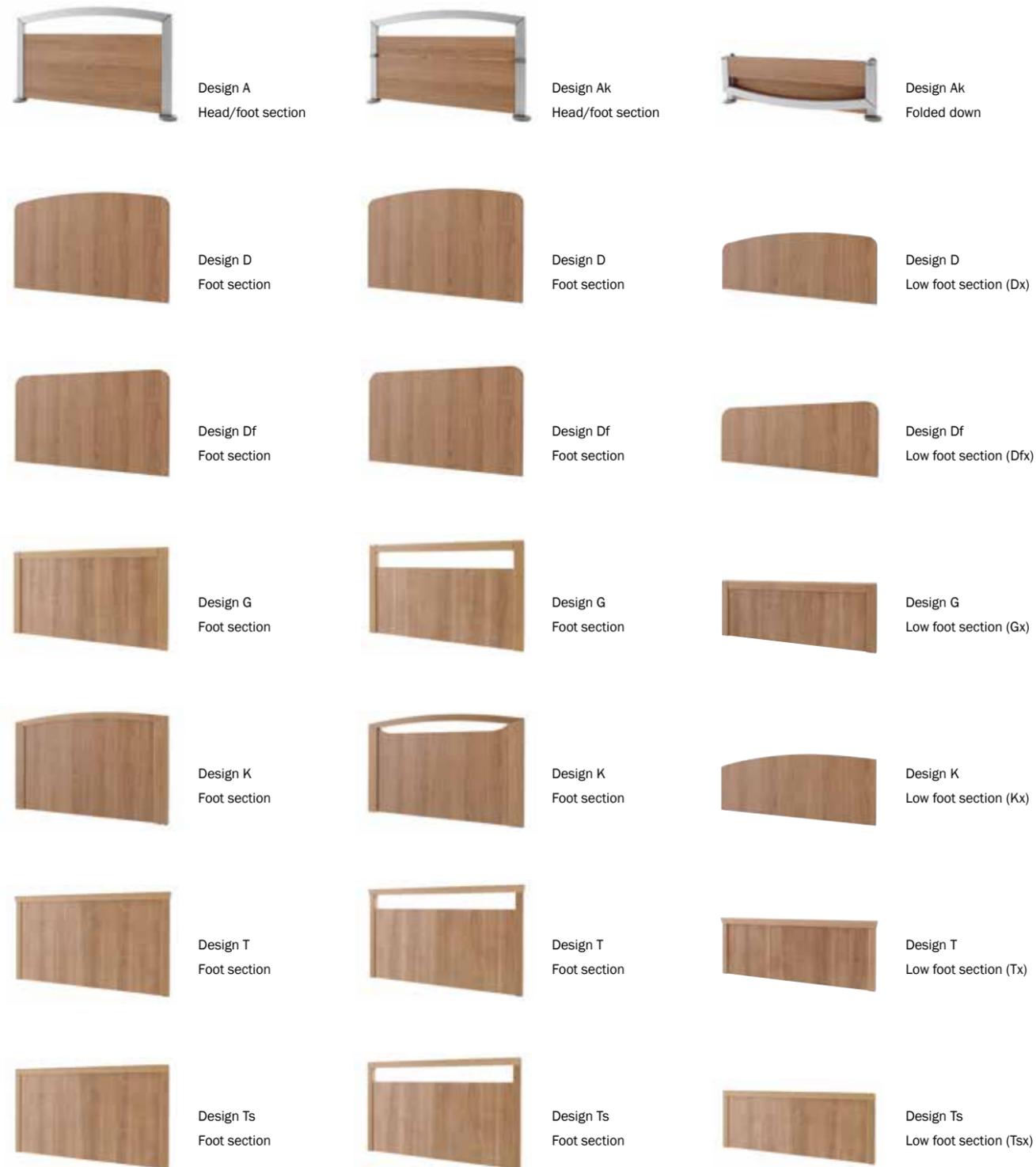

Colors | décors*




Color combination memoriana

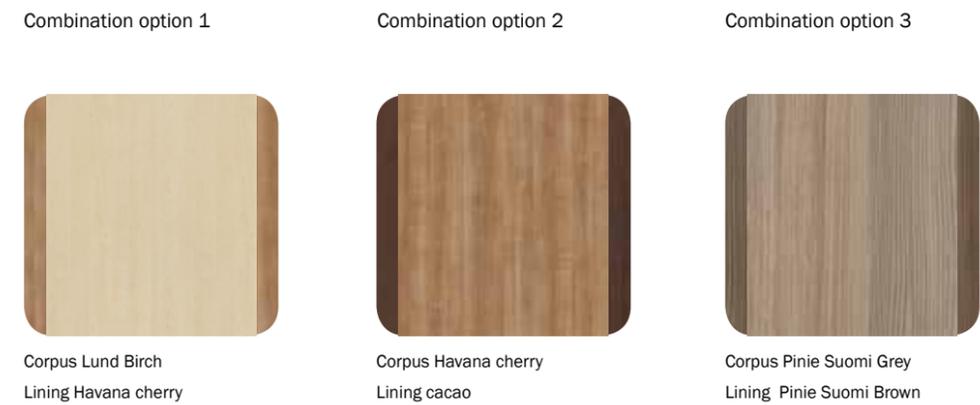


Furnishings and accessories

Dimensions

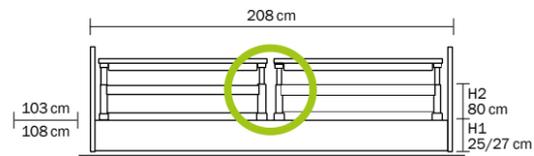


sentida with divided (10) side guards, without gap
Outer dimensions: L = 208 cm, W = 103/108 cm
Mattress dimensions: L = 200 cm, W = 87 cm
H1 25/27 cm (low position), H2 80 cm (highest position)

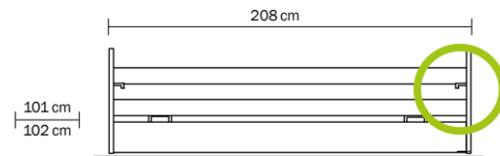


sentida with continuous 2-part
09 side guard with trigger release


Patient surfaces

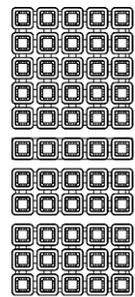

Comfort

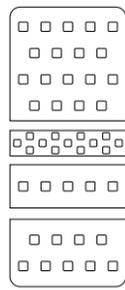


EasyClean

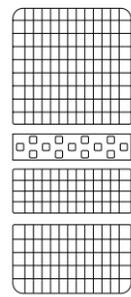


Standard | Aero


Bed extensions

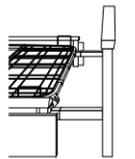


Temporary bed extension divided
10 side guard (standard)

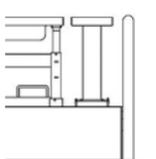


Permanent bed extension divided
10 side guard (optional)

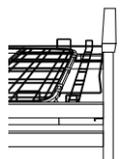

Permanent bed extension
09 side guard (optional)

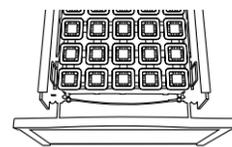

Bed extension with express release
(optional)

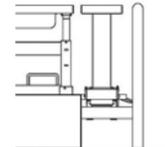

10-0570 / 10-0580
Protector for bed extension
by 20 cm (optional)

Options




Vertical wall stopper castors




3D horizontal and vertical wall stopper castors



CPR release on the backrest

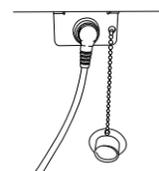


Plug & play use of the handset on one or both sides of the bed

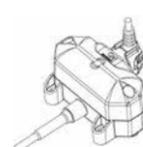


Underbed light (switchable)

Accessories

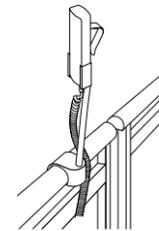


10-0419
Handset holder for divided side guards
10-0421
Handset holder, short

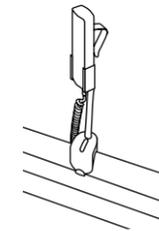


10-0422
Accessories and handset holder for continuous side guards
10-0421
Handset holder, short




10-0424
Accessories and handset holder for attachment to side panel
50-0164 urine bag holder (no picture)

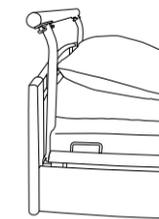

10-0534
Transport handle EasyMove

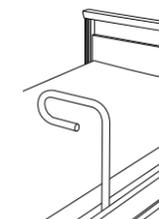


10-0715
Standing-up aid for 09 side guards as well as 3/4 side guards

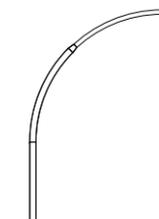


50-0237 soluna LED lamp with power supply unit
50-0550 soluna LED lamp without power supply unit
50-0547 soluna LED light with smooth lampshade without power supply unit




50-0206
Fall mat width 90 x length 200 cm, height approx. 4.5 cm

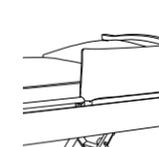


50-0174
Cushioning for divided 10 side guards

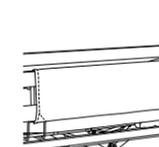


50-0338
Cushioning for 09 two-part side guards

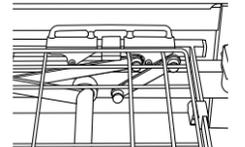


10-0563/10-0533/10-0521
Segufix holder
- for the thigh rest
- for the backrest
- for the seat section

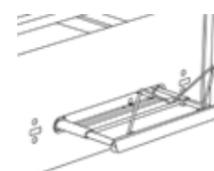


10-0617
Integrated design linen holder
Compatible with all beds (except 09 3-part side guards and A-kk design)

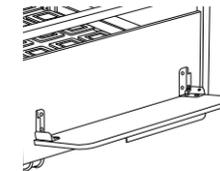

10-0719
Linen holder, compatible with all bed types

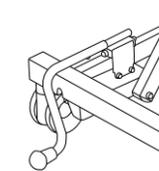


10-0707
Head end wall spacers, short (for sentida 5 + 6 only)

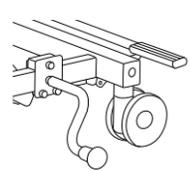

10-0708
Side-mounted wall spacer, short (for sentida 5 + 6 only)

Accessories – SafeSense® bed exit system

SafeSense® sensor mat 2

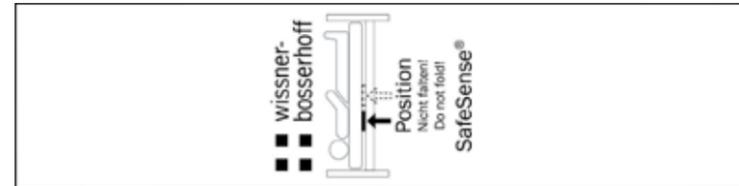


02-000876
SafeSense® sensor mat 2
80 x 20 cm

SafeSense® Bed-Exit-Box 2



04-1851
Digital bed exit box 2
cable or wireless version

SafeSense® connection kit

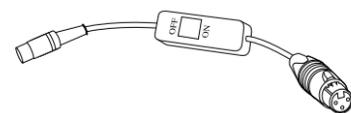

04-1508
Adapter cable with on/off switch and suitable plug for connection to the nurse call system (customer-specific) for the cable and wireless version



04-1608
Adapter cable without on/off switch with suitable plug to connect to the nurse call system (customer-specific) for the cable and wireless version

Our certificates – your safety

wissner-bosserhoff voluntarily has its production processes and management regularly audited by the TÜV. The company is certified in accordance with the regulations of the standard DIN EN ISO 9001:2008 for quality management systems as well as DIN EN ISO 13485:2003+AC:2007 for the manufacture of medical devices.



ISO 9001:2008 EN ISO 13485:2012 ISO 14001:2004 ISO 50001:2011



Inspection label for product certification and auditing



CB certificate Certificate for side guards

Design

wissner-bosserhoff places great emphasis on design. Not only do our products combine functionality and high-quality workmanship but also design at the highest level. They feature a high level of functionality and cozy design as well as convenient usability and ergonomics. Having won several prestigious design awards, sentida is now among the most important nursing home bed series on the international market.



Members of LINET Group

wissner-bosserhoff GmbH

Hauptstraße 4-6 | 58739 Wickede (Ruhr) | Germany

Phone: +49(0)2377.784-0 | Fax: -163 | info@wi-bo.de | www.wi-bo.de

